



## CONTAMINATION CONTROL

## Particle Counter Software



### Features and Benefits

- Enables 21 CFR Part 11 Compliance
- Configurable Security Privileges
- User Name and Password Access
- Lock Out After Failed Logins
- Password Aging
- Audit Trail
- Automatic Data Download
- Easy Data Export
- Data Integrity Checking
- Recipe Download
- Sample Location Name Download
- Create and Upload Recipes
- Create and Upload Sample Location Names

### TrakPro™ Lite Secure Download Software

TSI TRAKPRO Lite Secure download software for the range of AERO TRAK™ Portable Particle Counters enables compliance to FDA 21 CFR Part 11. Unique user names and passwords, audit trail, password aging, and lockout after failed logins, ensures secure, tamper proof data transfer. Individually assigned user security privileges enables managers to control levels of user access and activity.

TRAKPRO Lite Secure software automatically detects the AERO TRAK particle counter on connection to the computer and automatically downloads the data, location names and recipes, no user intervention is required. Data is then saved, printed, or exported. Exported data files can be verified for data integrity.

TRAKPRO Lite Secure software allows the user to securely modify, create, delete, and export AERO TRAK particle counter recipes and sample point names. Managing large quantities of sample point names and recipes when using multiple AERO TRAK particle counters could not be easier.



TRUST. SCIENCE. INNOVATION.

# Particle Counter Software

## Hardware (minimum requirements)

- Pentium 90 MHz, 32 MB RAM, 20 MB of hard drive space, or minimum requirements for the operating system, whichever is higher. Free hard drive space depends on the amount of records to be stored.
- A free USB port.

## Software (minimum requirements)

- Windows® XP operating system with Service Pack 2 or higher, Windows® Vista® (32-bit only) operating system with Service Pack 1 or higher, Windows® 7 operating system (32-bit or 64-bit).
- NET 2.0 or higher.


## Supported Instruments

- AEROtrak handheld particle counters: 9306 and 9303
- AEROtrak portable particle counters: 9310/9510/9350/9550 and 9500
- AEROtrak portable particle counter: 9110


## Ordering Information

To order TRAKPRO Lite Secure software, request TPL-S90.


Data Screen

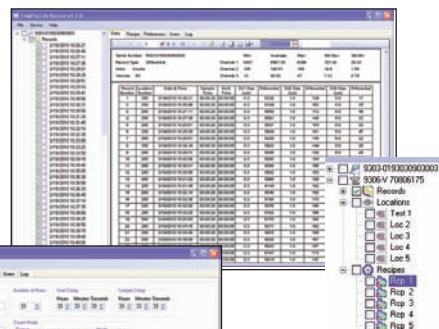

Recipe Screen

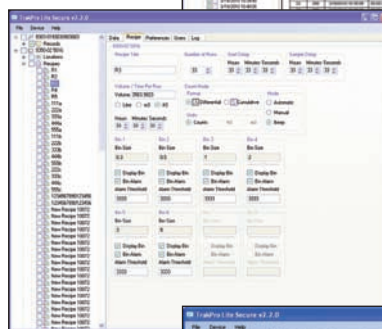

Log Screen

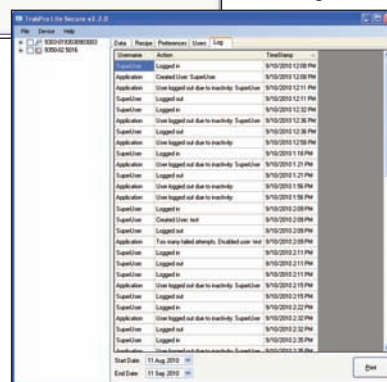

AeroTrak and TrakPro are trademarks of TSI Incorporated.  
Windows and Vista are registered trademarks of Microsoft Corporation.

Specifications are subject to change without notice.

**TSI Incorporated** - 500 Cardigan Road, Shoreview, MN 55126-3996 USA  
**USA** Tel: +1 800 874 2811 E-mail: info@tsi.com Website: www.tsi.com  
**UK** Tel: +44 149 4 459200 E-mail: tsiuk@tsi.com Website: www.tsiinc.co.uk  
**France** Tel: +33 491 11 87 64 E-mail: tsifrance@tsi.com Website: www.tsiinc.fr  
**Germany** Tel: +49 241 523030 E-mail: tsigmbh@tsi.com Website: www.tsiinc.de  
**India** Tel: +91 80 41132470 E-mail: tsi-india@tsi.com  
**China** Tel: +86 10 8251 6588 E-mail: tsibeijing@tsi.com  
**Singapore** Tel: +65 6595 6388 E-mail: tsi-singapore@tsi.com



TRUST. SCIENCE. INNOVATION.

Contact your local TSI Distributor or visit our website [www.tsi.com](http://www.tsi.com) for more detailed specifications.